


Gravimetric Certificate

Dr. Ehrenstorfer

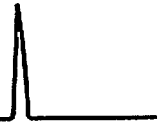

Product Identification

16659510 Pyriminobac-methyl (E)
Formula C₁₇H₁₉N₃O₆
Mol.Weight 361.34
CAS No. 147411-69-6

Reference Materials for Residue Analysis

Expiry Date 09.01.2017
Lot Number 21010CY
Store at 20°C in the dark

Please note: The expiry date is valid under recommended storage conditions only.

Gravimetric Data		Conc. (mg/l)	Purity %	Weight (mg)
Product Name Pyriminobac-methyl (E)		10.000	99.7	1.008
Solvent Information		Lot No.	Exact Quantity (ml)	
Solvent Cyclohexane		20801	100.00	
Traceability Data				
16659510 10119	neat product			
16659510 21010CY	10.000 mg/l			
Analytical Data		Method Details:		
Detection: HPLC/DAD	Acetonitrile:H ₂ O 9:1			
Column: ReproSil 100 C18 5µ 250x3				
Inj.-Vol.: 10.00 µl				
Flow: 1.0 ml/min				
Ret.-Time: 1.42 min.				
Identity check UV, RT				
Comment				
The uncertainty/tolerance of this standard is +/- 2.0 %, calculated in accordance with the EURACHEM/CITAC Guide - Quantifying Uncertainty in Analytical Measurement - Second Edition. The uncertainty given is the expanded combined uncertainty and represents an estimated standard deviation equal to the positive square root of the total variance of the uncertainty of components. The expanded uncertainty is U which is U _c (y)*K, where K is the coverage factor at the 95% confidence level (K=2). The expanded uncertainty is based on the combination of uncertainties associated with each individual operation involved in the preparation of this product.				

Certified on 09.01.2013
by C. Neukom



RES.NO.: 2274-01

The Laboratory Labor Dr. Ehrenstorfer-Schäfers is accredited by DGA as indicated by the Accreditation Certificate DGA-PL-4536.00 based on DIN EN ISO/IEC 17025:2005 for the weighing of amounts of substances for the preparation of standard solutions.

Labor Dr. Ehrenstorfer-Schäfers · Bgm.-Schlosser-Str. 6 A · 86199 Augsburg · Germany
Phone +49 821 906080 Fax +49 821 9060888 info@analytical-standards.com
The warranty for this product is limited to the purchasing price of this product.



DGA-PL-4536.00

=====


Acq. Operator : Dr. Heidrich
 Acq. Instrument : DAD1
 Injection Date : 10.10.2012 05:57:53
 Inj Volume : 10.000 µl

Seq. Line : 65
 Location : Vial 69
 Inj : 1

Acq. Method : D:\CHEM32\1\METHODS\91K.M
 Last changed : 06.09.2012 08:39:41 by Dr. Heidrich
 Analysis Method : V:\DAD1\METHODS\91K.M
 Last changed : 06.09.2012 08:39:41 by Dr. Heidrich
 Method Info : Acetonitrile : Water 9 : 1

Sample Info : Pyriminobac-methyl (E)

Additional Info : Peak(s) manually integrated
 DAD1A, Sig=205,2 Ref=off (V:\DAD1\201209KW41\16659510-8.D)

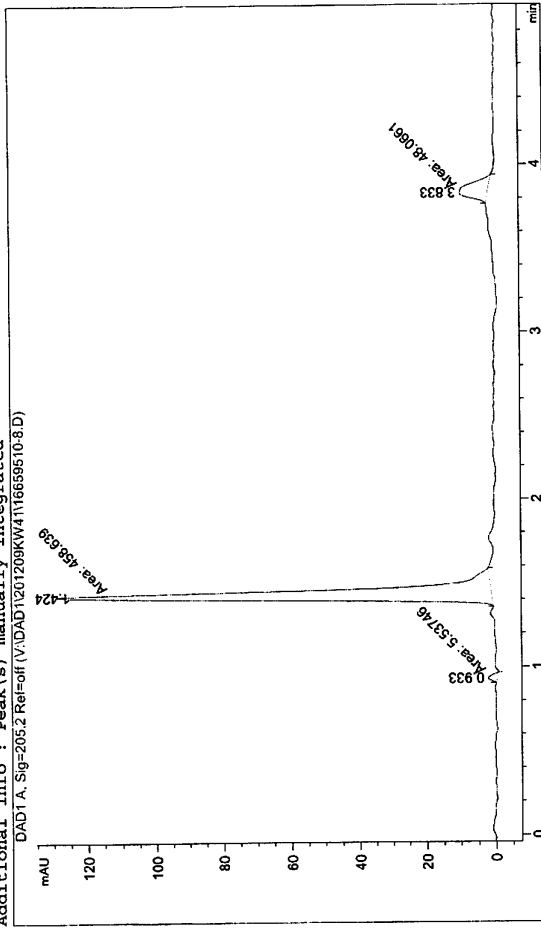

=====
 Area Percent Report
 =====

Sorted By	Retention Time	Area	Height	Area %
Multiplier:	: 1.0000			
Dilution:	: 1.0000			
Use Multiplier & Dilution Factor with ISTDs				

Signal 1: DAD1 A, Sig=205,2 Ref=off

Peak #	RetTime [min]	Sig	Type	Area [mAU*s]	Height [mAU]	Area %
1	0.933	1	MM T	5.53746	2.55851	1.0810
2	1.424	1	MM T	458.63882	128.83372	89.5355
3	3.833	1	MM T	48.06607	8.13473	9.3815

Totals : 512.24236 139.52696

Dr. Ehrenstorfer

Reference Materials for
Residue Analysis

1. Application:

This standard solution is designed for calibration or recalibration of chromatographic systems for the determination of the specified chemical compounds concerning identity and quantification. The product can also be used as reference material for interlaboratory studies to validate analytical procedures.

2. Raw material:

All raw materials used to prepare this standard solution are of the highest purity. After our production process each material is checked by several of the following methods, if applicable: UV-, IR-Spectroscopy, Elemental analysis, chromatographical properties (GC/FID, GC/ECD, HPLC/DAD, TLC) and physical properties (phase, colour, odor, melting point) and Karl-Fischer for detection of traces of water.

3. Manufacturing:

Prior to the production for each standard solution we calculate the necessary weight of the neat material, compensate the difference of the purity of the compound to 100%. We use an electronic scale capable of weighing to 0,000001 g with a built in automatic calibration function, which is executed minimum once a week. Once a month the balance is calibrated with weights complying with the OIML-IR-20 design requirements and traceable to the national prototype of the Physikalisch-Technische Bundesanstalt, Braunschweig, the German office of weights and measurements. Every year the balance is calibrated by the manufacturer service technician. This company is also certified by the DQS following ISO9001. All steps are documented conform to ISO 9001 requirements. The single components are weighed and filled into volumetric flask (class A glassware) with the exact quantity of solvent as indicated on the certificate. Solvents are stored and handled in special rooms which have the constant ambient temperature required from the manufacturer of the volumetric equipment. The tolerance of the weighing procedure and the dilution error adds to a maximum of +/- 1.0 %.

4. Packaging and Storage:

The final formulation was packaged in amber OPC ampoules to prevent photodegradation and then sealed by hand. Every reference material is stored under controlled condition. One sample of each lot is kept according GLP rules to bring our staff in a position to check the specific lot even years after the last item was sold.

5. Stability:

In regular intervals each lot is checked for stability. We guarantee the stability of the solution until expiry date given in the Gravimetric Certificate. In case, that the tests show a degradation within the period the customer will be notified by Dr. Ehrenstorfer. We recommend to store the ampoules in the dark at 20°C +/- 4°C

6. Gravimetric Certificate:

The documentation gives all the data of the production process with all the information necessary for traceability of each lot. Following GLP rules you are obliged to note the used product, lot number of the product, purity, exact weight/quantity of the product, name of the solvent, lot number of the solvent, exact quantity of the solvent, date of production, date of expiry and signature of the person in charge. Copies have to be authorized by stamp, signature of the person in charge according to your quality management handbook and date of authorization.

7. Analytical Quality Control:

The summary of the quality control procedure is documented in the separate Certificate of Analysis. Our company holds the Quality System Certificate DQS-Reg No.:2874-01 for the standard from the ISO 9001 / EN 29000 series and the scope as specified. The audit performed by the DQS has verified, that our quality system fulfills the requirements of DIN ISO 9001. The Laboratory Dr. Ehrenstorfer-Schäfers is accredited by DGA as indicated by the Accreditation Certificate DGA-PL-4536.00 based on DIN EN ISO/IEC 17025:2005 for the weighing of amounts of substances for the preparation of standard solutions.



The Laboratory Dr. Ehrenstorfer-Schäfers is accredited by DAR as indicated by the Accreditation Certificate DGA-PL-4536.00 based on DIN EN ISO/IEC 17025:2005 for the weighing of amounts of substances for the preparation of standard solutions.

Labor Dr. Ehrenstorfer-Schäfers · Bgm.-Schlosser-Str. 6 A · 86199 Augsburg · Germany
 Phone +49 821 906080 Fax +49 821 9060888 info@analytical-standards.com
 The warranty for this product is limited to the purchasing price of this product.



DGA-PL-4536.00