




Certificate of Analysis

Reference Standard

Pyridoxine Diacetate

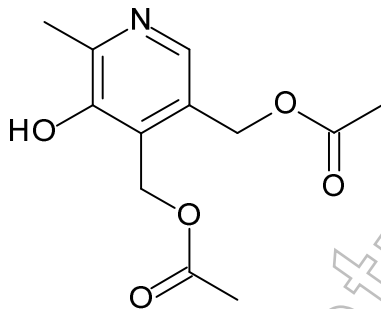

Molecular Formula: $C_{12}H_{15}NO_5$
Molecular Weight: 253.25
CAS Number: [26280-83-1]

Catalogue Number: MM0021.03
Lot Number: 42149
Long-term Storage: 2 to 8 °C, dark
Appearance: white solid
Melting Point: 107 °C
Assay 'as is': 99.4 %

Date of shipment: **2015-July-04**

This certificate is valid for two years from the date of shipment provided the substance is stored under the recommended conditions.




I. Identity

The identity of the reference substance was established by following analyses.

Ia. ¹H-NMR Spectrum

Conditions: 400 MHz, DMSO-d₆


The structure is confirmed with the signals of the spectrum and their interpretation.

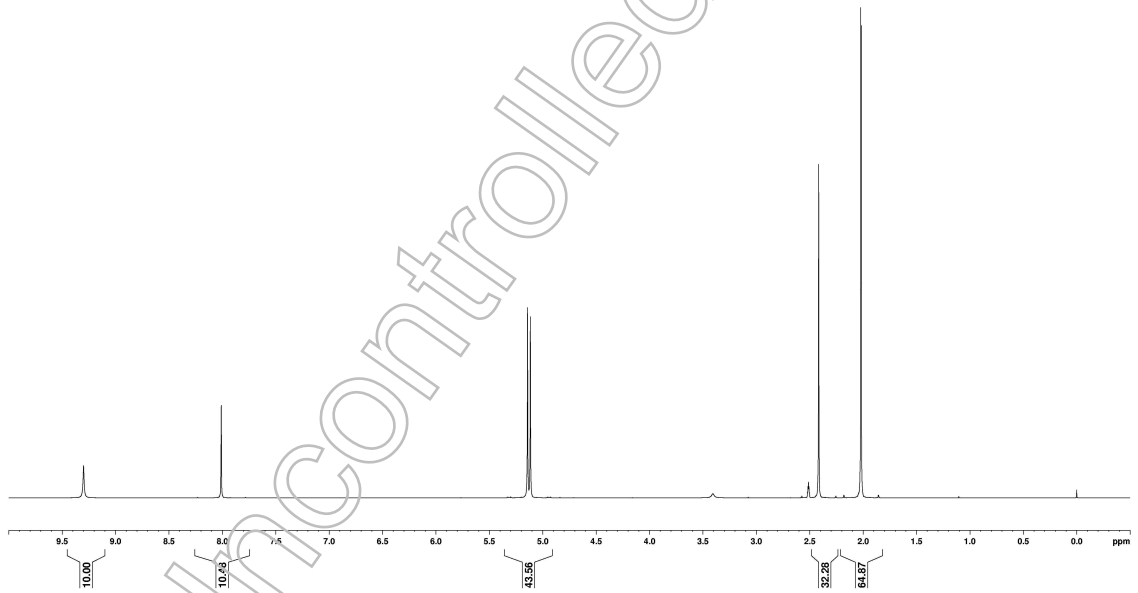



1b. Mass Spectrum

Method: 4.5 kV ESI; vaporization temperature: 200 °C, direct inlet

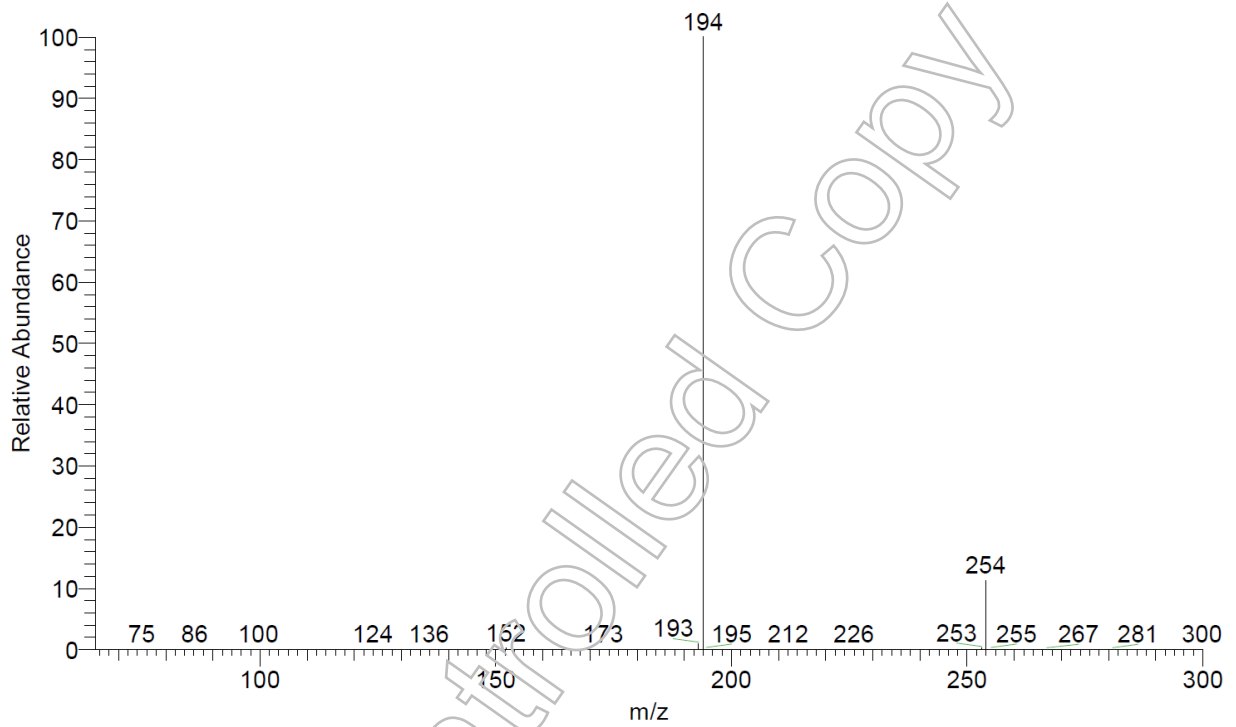


m/z	fragments
254	[MH]
194	[MH - CH ₃ COOH]

The signals of the mass spectrum and their interpretation are consistent with the structural formula.



Ic. IR Spectrum

Method: Attenuated Total Reflection Fourier Transform Infrared (ATR-FTIR) Spectroscopy

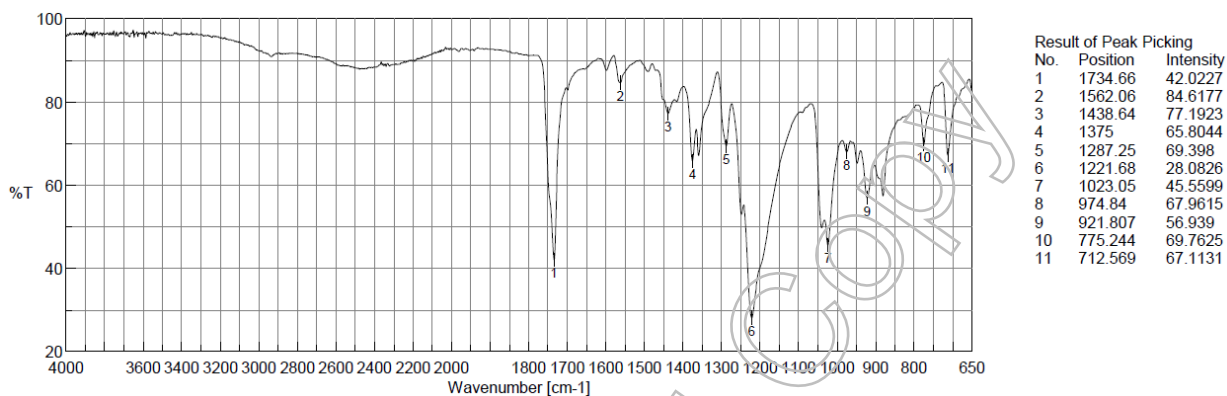

The signals of the IR spectrum and their interpretation are consistent with the structural formula.

II. Purity

The purity of the reference substance was analysed by high performance liquid chromatography (HPLC).

HPLC Conditions:

Column:

RP 60 Select B
5 μ m, 125 x 4 mm

Conditions:


1.0 ml/min, 40 °C
0 – 5 min Water/Acetonitrile 90/10
5 – 8 min Water/Acetonitrile to 50/50
8 – 11 min Water/Acetonitrile to 90/10
11 – 16 min Water/Acetonitrile 90/10 (v/v);
0.1 % H₃PO₄

Detector:

DAD
295 nm

Injector:

Auto
1.5 μ l; 0.316 mg/ml in
Acetonitrile

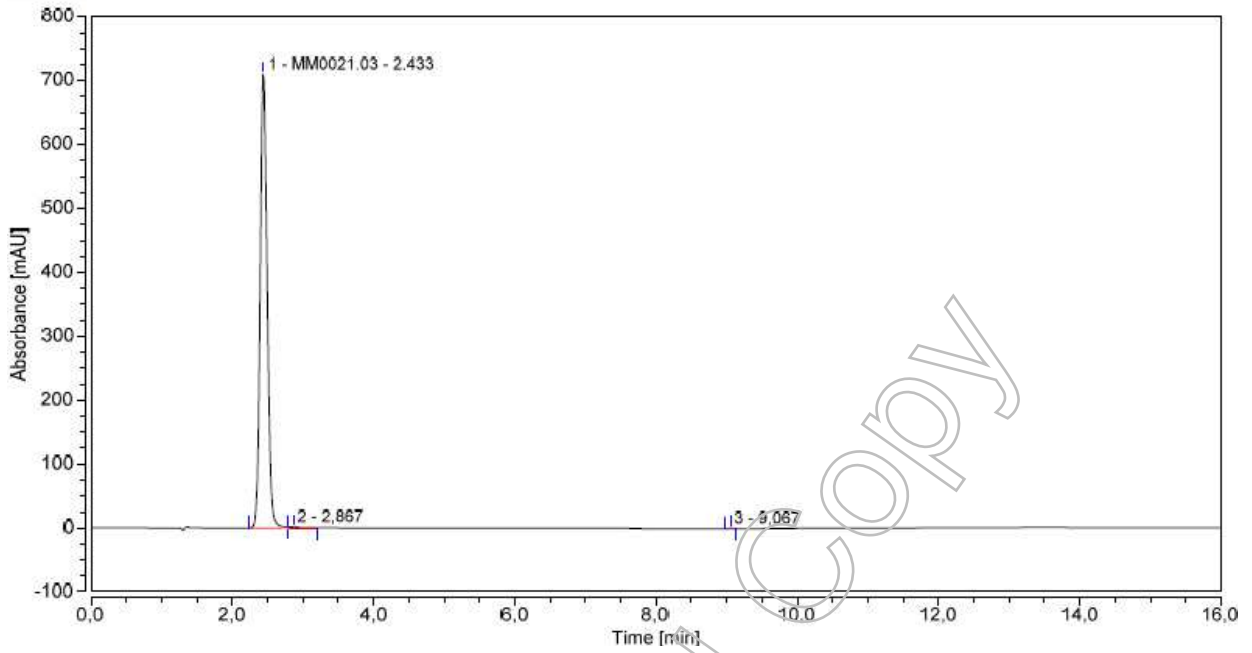

Area Percent Report - Sorted by Signal

Pk #	Retention Time	Area	Area %
1	2.433	82.033	99.49
2	2.867	0.329	0.40
3	9.067	0.091	0.11
Totals		82.453	100.00

For the calculation the system peaks were ignored. The content of the analyte was determined as ratio of the peak area of the analyte and the cumulative areas of the purities, added up to 100 %.

Results:

Average 99.49 %
Number of results n=3
Standard deviation 0.01 %



III. Water Content

Method: Karl Fischer titration

Results:

Average	0.10 %
Number of results	n=3
Standard deviation	< 0.01 %

IV. Residual Solvents

Method: ¹H-NMR

No significant amounts of residual solvents were detected (< 0.05 %).

V. Final Result

Total impurities (HPLC)	0.51 %
Water content	0.10 %
Residual solvents	n. d. (not detected)
Assay (100 % method) ¹	99.39 %

The assay is assessed to be 99.4 % 'as is'

The assay 'as is' is equivalent to the assay based on the not anhydrous and not dried substance respectively.

Release Date: 2013-11-27

LGC GmbH

Dr. Sabine Schröder
Product Release

¹ The calculation of the 100 % method follows the formula:

$$\text{Assay (\%)} = (100 \% - \text{KF} - \text{RES}) \times \frac{\text{Purity HPLC (\%)}}{100 \%}$$

Water (KF) and Residual solvents (RES) are considered as absolute contributions, HPLC purity is considered as relative contribution.