

15-ACETYLDEOXYNIVALENOL IN ACETONITRILE

1. General information

This document is designed and the certified value(s) and uncertainty(ies) are determined in accordance with ISO Guide 31 [1], ISO Guide 35 [2] and Eurachem / CITAC Guides [3,4].

2. Description of the Reference Material (RM)

Name:	15-Acetyldeoxynivalenol in Acetonitrile
CAS number:	88337-96-6
Catalog number:	B-MYC0360-5, B-MYC0360-1, B-MYC0360-2
Lot #:	L13374A
Expiry date:	13.03.2015
Starting material:	15-Acetyldeoxynivalenol, Lot # S12411A, Romer Labs Diagnostic GmbH
Physical description of RM:	Solution of 15-Acetyldeoxynivalenol in Acetonitrile
Packaging and amount of RM:	<u>B-MYC0360-5, B-MYC0360-2</u> : Amber glass ampoules fitted with teflon faced butyl septa and PP screw caps, solution of 5 mL and 2 mL <u>B-MYC0360-1</u> : Amber glass ampoules fitted with teflon faced butyl septa and aluminium crimp cap, solution of 1 mL
Name and address of the manufacturer:	Romer Labs Diagnostic GmbH Technopark 1, 3430 Tulln, Austria www.romerlabs.com
Name and address of the supplier:	LGC Standards GmbH Mercatorstraße 51, 46485 Wesel, Germany Tel +49(0)2 81 98 87 0, Fax +49(0)2 81/98 87 199 www.lgcstandards.com

2.1 Intended use of the RM

- for laboratory use only
- calibration of analytical instruments

2.2 Instruction for the correct use of the RM

The ampoules should be stored at 2-8°C or below in a dark place. Before usage of the RM, the ampoules should be allowed to warm to room temperature. The recommended minimum sub-sample amount for all kinds of application is 100 µL. The expiry date of this RM is based on the current knowledge and holds only for proper storage conditions in the originally closed flasks/packages.

2.3 Hazardous situation

The normal laboratory safety precautions should be observed when working with this RM. Further details for the handling of this RM are available as safety data sheet (SDS).

Hazardous Ingredients	Concentration in %	Pictograms	Signal word	Hazard statement(s)
Acetonitrile	> 99.9		Danger	H225, H302, H312, H319, H332

3. Certified values and their uncertainties

15-Acetyldeoxynivalenol in Acetonitrile		
Compound	Mass concentration ^a	
	Certified value ^b	Uncertainty ^c
15-Acetyldeoxynivalenol	101.0 µg/mL	± 0.6 µg/mL

^a Values are based on preparation data and confirmed experimentally by HPLC-UV
^b Mass concentration based on weighed amount, purity and dilution step
^c Expanded uncertainty U (k = 2) of the value u_c according to GUM [5]

3.1 Calculation of uncertainty

The uncertainty of the calibrant solution was calculated on the basis of preparation [6].

Uncertainty components	Description	Standard uncertainty (u)	
Purity (P) of solid 15-Acetyldeoxynivalenol	P = 99.7 % ± 0.3 %	u (P) = 0.2 %	a
Weighing procedure weighted sample: m _{ws} = 25.332 mg	repeatability: 0.03 mg linearity: 0.012 mg	u (m) = 0.03 mg	b
Dilution procedure volumetric flask: V _f = 250 mL	calibration: 250 mL ± 0.15 mL repeatability: 0.03 mL volume expansion solvent	u (cal) = 0.06 mL u (rep) = 0.03 mL u (Vol. exp.) = 0.6 mL u (V) = 0.6 mL	c d e f

^a Maximum tolerance of purity (rectangular distribution) was divided by $\sqrt{3}$

^b Estimation of this u-value is based upon the values for repeatability and linearity described in the user manual of the microbalance

^c A triangular distribution (division by $\sqrt{6}$) was chosen for the calculation of u (cal)

^d Based on a series of ten fill and weigh experiments on a typical 250 mL flask; the value was used directly as a standard deviation

^e Based on the density of 0.7857 g/cm³ at temperature T = 20°C and a maximum temperature variation of ± 3°C, of volume expansion, relative volume expansion coefficient of acetonitrile is $1370 \cdot 10^{-6}/^{\circ}\text{C}$ [7], volume expansion term (rectangular distribution) was divided by $\sqrt{3}$

^f The three contributions are combined to give the u (V) = $\sqrt{u(\text{cal})^2 + u(\text{rep})^2 + u(\text{Vol. exp.})^2}$

Calculation of the combined uncertainty u_c and the expanded standard uncertainty U for 15-Acetyldeoxynivalenol

$$c_{\text{Toxin}} = \frac{10 \times m_{\text{ws}} \times P}{V_f} = \frac{10 \times 25.332 \times 99.7}{250} = 101.0 \text{ mg / L}$$

$$\frac{u_c(c_{\text{Toxin}})}{c_{\text{Toxin}}} = \sqrt{\left[\frac{u(P)}{P}\right]^2 + \left[\frac{u(m)}{m_{\text{ws}}}\right]^2 + \left[\frac{u(V)}{V_f}\right]^2} = \sqrt{\left[\frac{0.2}{99.7}\right]^2 + \left[\frac{0.03}{25.332}\right]^2 + \left[\frac{0.6}{250}\right]^2} = 0.003$$

$$u_c(c_{\text{Toxin}}) = c_{\text{Toxin}} \times 0.003 = 101.0 \times 0.003 = 0.3 \text{ mg / L}$$

Calculation of expanded standard uncertainty U using a coverage factor k = 2

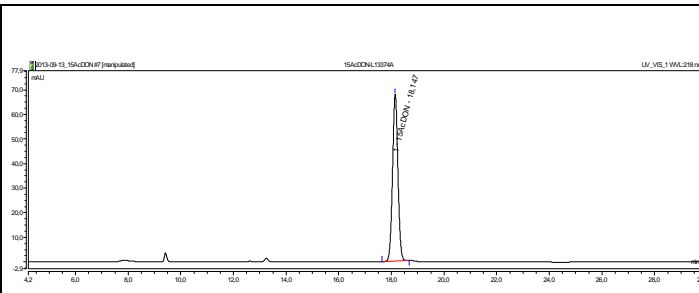
$$U(c_{\text{Toxin}}) = u_c(c_{\text{Toxin}}) \times 2 = 0.3 \times 2 = 0.6 \text{ mg / L} = 0.6 \mu\text{g / mL}$$

4. Discussion of traceability

This calibrant is certified on the basis of gravimetric preparation [6]. Thus the certified value (mass concentration of 15-Acetyldeoxynivalenol) is based on the weighed amount of the starting material and is therefore traceable to the stated purity of the solid mycotoxin. High purity material represents a practical realization of concentration units, through conversion of mass to molar quantity.

5. Confirmation of certified value by HPLC-UV

The certified concentration of 15-Acetyldeoxynivalenol of the gravimetric prepared solution was confirmed by HPLC-UV against an independently prepared reference batch of 15-Acetyldeoxynivalenol calibrant Lot# L13092D.

column	Phenomenex Luna C18(2), 250x3.00 mm, 5µ		 <p>Figure 1: HPLC-UV chromatogram of 15AcDON calibrant Lot # L13374A</p>								
injection volume	100 µL sample										
solvent A	Acetonitrile / Water (10 / 90)										
solvent B	Acetonitrile										
sample dilution	1:20 with water										
flow rate	0.5 mL / min										
gradient	time in minutes (min)	% solvent B									
	0 – 3	0									
	3 – 6	0 - 11									
	6 – 20	11									
	20 – 21	11 - 0									
	21 – 30	0									
DAD settings	218 nm										
			<table border="1"> <thead> <tr> <th></th> <th>time [min]</th> <th>area</th> <th>concentration ^a [µg/mL]</th> </tr> </thead> <tbody> <tr> <td>15-Acetyldeoxynivalenol</td> <td>18.147</td> <td>16.049</td> <td>101.2 ± 2.5</td> </tr> </tbody> </table>		time [min]	area	concentration ^a [µg/mL]	15-Acetyldeoxynivalenol	18.147	16.049	101.2 ± 2.5
	time [min]	area	concentration ^a [µg/mL]								
15-Acetyldeoxynivalenol	18.147	16.049	101.2 ± 2.5								
			^a Mean of 7 replicate measurements against reference batch, confidence interval with P = 95 %								

6. Further information

The purchaser must determine the suitability of this product for its particular use. LGC Standards GmbH makes no warranty of any kind, express or implied, other than its products meet all quality control standards set by LGC Standards GmbH. We do not guarantee that the product can be used for a special application.

approved for release by: Dipl. Chem. Ulrich Berger

date: 13.09.2013

This document has been computer generated and is valid without a signature.

References:

- [1] ISO Guide 31, 1-7, (2000), "Reference materials – Contents of certificates and labels"
- [2] ISO Guide 35, 1-7, (2000), "Certification of Reference Materials – General and Statistical Principles"
- [3] Eurachem / CITAC Guide, 1-37, (2003), "Traceability in Chemical Measurement"
- [4] Eurachem / CITAC Guide, 1-120, (2000), "Quantifying Uncertainty in Analytical Measurement"
- [5] International Organization for Standardization (ISO), (1995), "Guide to the expression of uncertainty in measurement", 1st Ed. Geneva, Switzerland
- [6] R.D. Josephs, R. Krska, S. MacDonald, P. Wilson, H. Pettersson, J. AOAC Int. **86**, 50-60, (2003), "Preparation of a Calibrant as Certified Reference Material for Determination of the Fusarium Mycotoxin Zearalenone"
- [7] E.W. Flick, (1996), "Industrial Solvents Handbook", 3rd Ed., Noyes Data Corp. Westwood NJ